# RotoFlex™ Plus Tube Rotator

Operator's Manual

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
Manuale d'uso

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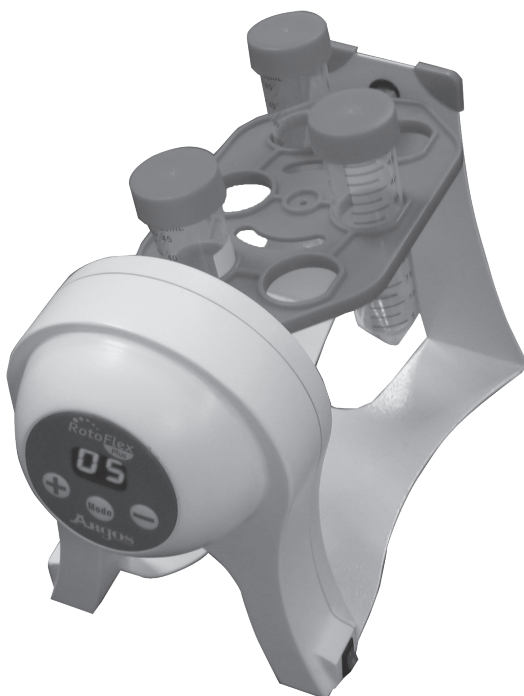**ARGOS**  
TECHNOLOGIES *manufacturing solutions for science*

**Argos Technologies, Inc.**

[customer-service@argos-tech.com](mailto:customer-service@argos-tech.com)

Toll Free: 800-886-8675

International: 847-622-0456

[www.argos-tech.com](http://www.argos-tech.com)

## Introduction

The RotoFlex Plus has been designed for applications that require a more rigorous mixing action than our standard RotoFlex instrument offers. In addition to a faster rotational speed (adjustable between 10 -40 rpm) versus 22 rpm on the standard RotoFlex; the RotoFlex Plus can also be operated in oscillating mode for a more complete mixing action. During the oscillating cycle, the rotisserie rapidly vibrates back and forth 9 times before rotating 360°. This reciprocating motion ensures a thorough mixture of all tubes; similar to a mild vortexing action.

To change rotisserie, simply snap on the desired tube size to the support bar. No tools required.

---

## The RotoFlex Plus is ideal for

- Hybridizations
  - DNA purification
  - Sample mixing / resuspensions
  - Protein binding
  - Immunoprecipitation
  - Any procedure where thorough mixing is required
- 

## Package Contents

Tube Revolver, Cat. Number <b>R2200</b> .....	1 Each
Rotisserie for 0.5 - 1mL tubes, Cat. Number <b>R2002</b> .....	2 Each
Rotisserie for 1.5 - 2mL tubes, Cat. Number <b>R2003</b> .....	2 Each
Rotisserie for 5 - 7mL tubes, Cat. Number <b>R2004</b> .....	2 Each
Rotisserie for 10 - 15mL tubes, Cat. Number <b>R2005</b> .....	2 Each
Rotisserie for 50mL tubes, Cat. Number <b>R2006</b> .....	1 Each
Rotisserie Axle, Cat. Number <b>R2007</b> .....	1 Each
Rotisserie Axle, for 50mL tubes, Cat. Number <b>R2008</b> .....	1 Each
O-Ring Pack (10 each), Cat. Number <b>R2009</b> .....	1 Pack
Replacement Fuse, Cat. Number <b>R2010</b> .....	1 Each
Power Supply, Cat. Number <b>R2011</b> .....	1 Each
Operator's Manual .....	1 Each

---


## Specifications

<b>Model</b>	R2200
<b>Rotation Speed</b>	10-40 rpm
<b>Capacity</b>	42 x 0.5-0.8mL tubes / 30 x 1.5-2.0mL tubes / 14 x 5-7mL tubes / 14 x 15mL tubes / 6 x 50mL tubes
<b>Ambient Temperature</b>	4°C to 60°C
<b>Speed Control</b>	Variable
<b>W x D x H (cm)</b>	26 x 14.8 x 19.5
<b>Power</b>	AC 100V – 240V, 50Hz / 60Hz
<b>Noise Level</b>	Noiseless Operation
<b>Weight</b>	2.4 lb / 1.09kg
<b>Certifications</b>	CE Marked

**This device complies with RoHS 2009/95/EC**

## Panel Instruction

1. Increase Speed
2. Decrease Speed
3. Selection of operation mode
4. LED to display settings and working status
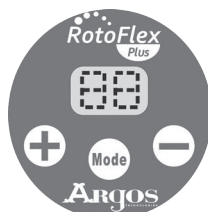
## Operating Instructions

1. Make sure the power switch is turned off before connecting the plug to a grounded wall socket.
2. Place the RotoFlex Plus on a level, horizontal and secure surface.
3. Select the rotisserie required for your application (see page 4) for rotisserie installation.
4. Insert power supply pin into the jack located on the base of the leg that houses the digital display. Insert the UL approved power cord into the grounded wall socket using the proper adapter supplied with the RotoFlex Plus.
5. Load the tubes and adjust the rotisserie to the desired angle. Be sure to distribute tubes evenly to achieve a smooth rotating motion.
6. Start the rotator by pushing down the power switch located on the leg opposite the power supply jack.
7. Rotational Mode:
  - a. Push the “mode” key to select the rotational mode. The digital display will be activated and register 10 rpm initially.
  - b. By using the up and down arrow keys located below the digital display, you can adjust the speed of the RotoFlex between 10-40 rpm.
  - c. The digital LED display will show the last speed setting when in rotational mode.

## 8. Oscillating Mode:

- Push the “mode” key to select the oscillating mode. The digital display will be activated and register 00.
- The rotisserie will oscillate back and forth 9 times in a cycle while rotating at a rate of 3 rpm. Once the oscillation cycle is complete, the rotisserie will begin to rotate at a speed of 40 rpm until the next oscillation cycle occurs, every 2 seconds or one complete rotisserie revolution.
- The up and down arrow keys will not function in the oscillation mode.
- By depressing the “mode” key once again, the instrument returns to the rotational mode, displaying the last speed selected.


## 9. When finished, turn the power switch off before removing tubes.


**Caution:** Do not touch plug with wet hands or use the device near a water source.

## Changing Rotisseries

- Firmly grasp the RotoFlex Plus rotator with your left hand while holding the rotisserie axle with your right hand. Gently pull the rotator body with your left hand outward until the axle is free of the unit.

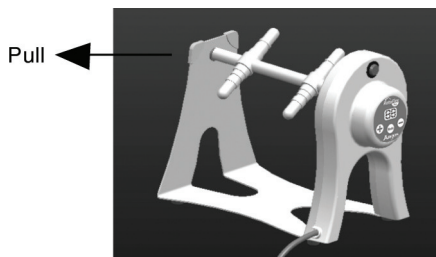- Remove the opposite side of the axle body by pulling it gently from the RotoFlex Plus housing.
- Select the rotisserie you wish to utilize and mount them on the RotoFlex Plus axle.

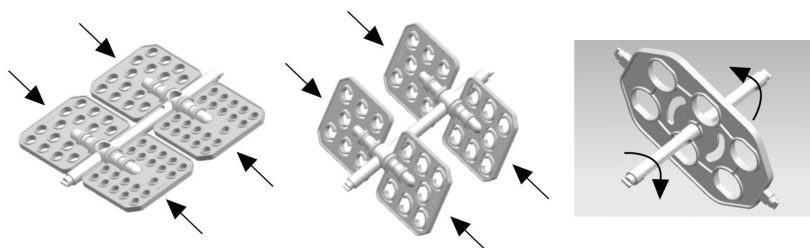- Insert the axle back into the RotoFlex Plus housing and turn until the axle clicks into place.
- Gently pull the rotator body (as in step 1) and reinsert the axle.
- Turn the unit on to allow the rotisserie axle to align with the RotoFlex Plus motor.

## Maintenance & Cleaning

- No lubrication is required.
- Before cleaning be sure to unplug the power cord from the wall socket.
- Remove the rotisseries.
- The housing can be cleaned using a soft damp cloth and a mild non-corrosive detergent.
- After cleaning, ensure that all parts are thoroughly dry before attempting to operate the unit.
- **Do not immerse the RotoFlex Plus in liquid.**

### Warnings

- Not to be used to mix flammable samples.
- Avoid aggressive solvents or vapors.
- Take appropriate precautions when mixing hazardous samples.

## Troubleshooting

Error	Cause	Solution
Unit does not start	Power cable not connected	Check that power cable is firmly plugged into wall outlet.
	Power Supply	Check that there is power to the wall outlet.
	Power Switch	Check that power switch is at ON status.
If suggested solutions are unsuccessful, please contact Argos Technologies Technical Support at Argos Technologies 800-886-8675 (Inside U.S.) 847-622-0454 (Outside U.S.) <a href="http://www.argos-tech.com">www.argos-tech.com</a>		

## Warranty

Argos RotoFlex Plus Tube Rotator is warranted to be free of defect in materials and workmanship for a period of 12 months from the date of purchase, providing that:

- The unit has not been opened or attempted to be opened by the end user.
- The user ensures that it has been carefully cleaned and decontaminated.
- The unit has been operated according to the instructions and not in any improper way.

In order to qualify for warranty repair; the customer must return the warranty card supplied with the unit that includes all of the information required to validate proof of purchase. Should warranty service be required, repair or replacement will be provided by Argos without charge. Freight shipping the unit to Argos will be the responsibility of the purchaser. There are no other warranties or guarantees, implied or otherwise expressed including any warranties of application or merchantability.

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Copyright © 2012 Argos Technologies, Inc.